ELECTRIC
HEATING PAD

GM-90

EN ELECTRIC HEATING PAD GM-90 USER'S MANUAL

- Two temperature levels
- Soft fabric cover removable
- Double safety thermostat


IMPORTANT INSTRUCTION FOR END-USERS:

1. Read the following instruction carefully before using the heating pad and keep it to ensure a proper future use.
2. This pad has a switch with a sliding key that admits two well-defined positions. These positions allow a temperature selection adaptable to your needs. In the second position of the key, the temperature is maximal. It is advised to start using the pad in the second position, and once the pad starts to heat, select the first or leave the second position (depending on your need).
3. Before each use carefully check the inner cover and the wire to make sure there are no damages or worn out pieces. If they are detected, do not use the heating pad and contact the after sales service center.
4. In order to get an excellent long-term performance and to avoid risks arising from the product misuse, we highly recommended to follow these instructions for use:
 - Do not use the heating pad if it is wet. Bear in mind it operates on household current.
 - Do not sit on the pad, do not crush or fold it, (place the pad where heat is applied on the body and not under it).
 - Do not use it without the cloth cover, do not pull the pad by the cord and do not use the cord as a handle.
 - Never use pins, needles or any other sharp metal object to hold the pad into the desired position.
 - Do not use it on children, handicapped or sleeping people or those with poor blood flow.
 - Do not use the pad while sleeping.
5. To avoid burning please follow the next instructions:
 - Do not cover it with blankets (burns may be caused regardless of control setting).
 - Do not apply cream or liniment; they may cause more heat than the recommended.

Do not use the pad at the maximum temperature settings. Check the skin under the pad frequently to avoid burns and blisters. Unplug the pad after being used – without exception. The device must not be used by people who are mentally/physically unable, with reduced sensory capacity, children unless they are supervised or instructed by a person in charge of their safety. Children must be supervised while using device in order to prevent them from playing with it.

CLEANING AND MAINTENANCE:

Remember to remove the cover from the pad before this operation. Put the heating pad away without pressing on it and avoiding crumples. Dry-cleaning is allowed.

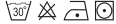TECHNICAL SPECIFICATIONS:

- 220V-240V/~50-60Hz/45W
- Size: 30 x 40 cm
- Weight: 520 g
- Composition: inner cover - 100% PVC, cloth cover - polar fleece.

■ Shijiazhuang Beijiren Electric Appliance Co., Ltd. South of Yishahe Bridge Chenganzhen, Xinle City, Hebei Province, China

Technical alternations reserved.

Under condition of the strict following to the operation instruction: warranty - 1 year from date of purchase. Operating life - 4 years.

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